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Please send a money order in the amount of \$9,95 (for return postage) should you need to have your scale repaired.

TOTAL BODY[®]

MASDEL COLLECTION^{MC/TM}

Dear Customer:

Thank you for purchasing a TOTAL BODY electronic bathroom scale. We hope you will find it satisfactory. Please read these instructions carefully to take full advantage of your scale.

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
■ BATTERY

Type. 1x3V lithium battery (CR2032).
Please remove the isolation sheet from the battery contact before use.

■ TO REPLACE BATTERY

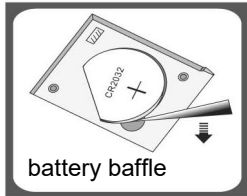Picture

Condition : Open battery case cover on the back of scale, if battery case is the type shown in picture,

- Take out the used battery with the aid of a sharp object, as shown in picture.
- Install new battery by putting one side of the battery beneath the battery baffle and then pressing down the other side.

■ OPERATION


STEP 1. Put the scale on a hard & flat surface (avoid carpet or soft surface).

STEP 2. Gently step on the scale, then the scale will be automatically turned on. Stand evenly until on the scale without moving and wait until your weight shown on display is stable and locked.


■ AUTOMATIC SWITCH OFF

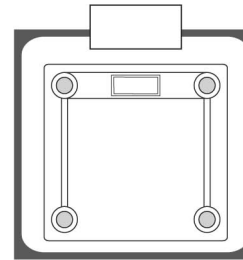When you step off the scale, it will automatically switch off. Auto shut-off occurs if display shows "0.0" or show the same weight reading for approximate 8 seconds.

■ TO CONVERT WEIGHT UNIT

If the scale is with both measuring system of Metric (kilograms) and Imperial (pounds), you can convert weight unit as following:

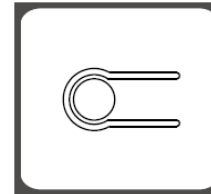Condition : if the unit button is on the back of the scale,

- Find unit button on the back of scale.
- Press the unit button to convert scale unit between kg-lb.

■ OVERLOAD INDICATION

Scale overload. (max. weight capacity is indicated on the scale).
Please step off to avoid damage.

■ LOW BATTERY INDICATION

Battery power is low. Please install new battery.

■ RE-TEST

Display " L " when you step off the scale, means this weight reading with error. Please re-test.

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ADVISE FOR USE & CARE

1. The glass platform will be slippery when it is wet. Keep it dry!
2. Keep standing still during the course of weighing.
3. DO NOT strike, shock or drop the scale.
4. Treat the scale with care as it is a precise instrument.
5. Clean the scale with a damp cloth and avoid water getting inside the scale.
DO NOT use chemical/abrasive cleaning agents.
6. Keep the scale in cool place with dry air.
7. Always keep your scale in horizontal position.
8. If failed to turn on the scale, check whether battery is installed, or battery power is low. Insert or replace new battery.
9. If there is error on display or cannot turn off the scale for a long time, please take off the battery for about 3 seconds, then install it again to remove soft-fault.
If you cannot resolve the problem, please contact your dealer for consultation.
10. Do not use it for commercial purpose.
11. **1 year warranty.**